




Yes, you can.™

The Invacare® MNS400 is a full-depth (8") powered alternating pressure air mattress replacement system which reduces interface pressure for the treatment of Stage I to Stage IV pressure ulcers. The MNS400 reduces interface pressures below capillary closure by alternating pressure in the 8-inch air cell bladders every 5, 10 or 15 minutes. The power unit control combines simplified operations with effective results.



### FEATURES

- Portable, lightweight mattress replacement system fits most standard medical beds
- Single pump 50 liter per minute air flow
- 16 Individual 8" nylon air cells reduce interface pressure below capillary closure while providing low interface pressure over the entire alternating cycle
- Mattress cells can be set to alternate every 5, 10 or 15 minutes to reduce interface pressures
- Waterproof, anti-microbial Ventex Recovery 5 therapy cover reduces friction and shear
- Durable ballistic nylon bottom protects the mattress cell from any abrasive surface
- Lower safety mattress supports the patient in the event of power failure
- Durable lower transport mat for patient transport or power outages

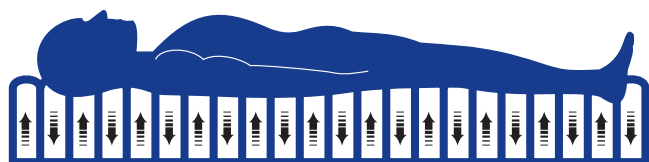

### Alternating Mode

In alternating mode, the 16 adjacent air cells of the Invacare MNS400 alternately inflate and deflate at 5, 10, or 15 minute adjustable intervals, which redistributes the pressure against the skin to promote capillary circulation.



### User-Friendly Controls

- 10 pressure settings for patient comfort range
- Power failure alarm
- Power unit is lightweight and compact
- Soft touch, easy wipe, sealed keyboard
- Low pressure alarm
- Dual connector rapid attach and CPR deflate at the power unit

### Specifications

Model no. MNS400-S (Standard width)  
Model no. MNS400-N (Narrow width)  
Model no. MNS400-B (Bariatric width)

Electrical power	115vac/60Hz
Dimensions	
Power unit	11" W x 12.5" H x 5.25" D
Standard Mattress	80" L x 36" W x 8.5" H
Narrow Mattress	80" L x 32" W x 8.5" H
Bariatric Mattress	80" L x 42" W x 8.5" H
Weight	
Mattress	
Standard mattress	22 lb.
Narrow mattress	20 lb.
Bariatric mattress	24 lb.
Power unit	9 lb.
Weight capacity	
MNS400-S	350 lb.
MNS400-N	350 lb.
MNS400-B	600 lb.
Material	
Mattress cells	70 denier nylon with a polyurethane coating CA# 117 Fire Retardency
Therapy pad	70 denier nylon top surface with a polyester backing for comfort CA# 117 Fire Retardency
Safety code approval	CE, UL2601, CSA, ETL
Limited warranty	
Power unit	1 year
Mattress	6 months

Invacare Corporation  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(866) 465-4387

Invacare and the blue medallion design are registered trademarks and Yes, you can. is a trademark of Invacare Corporation. Specifications are subject to change without notification.

Manufactured for Invacare Corporation by Genadyne Biotechnologies, Inc.

© 2005  
Invacare Corporation

Form No. 04-216



"Invacare says 'Yes, you can.' and I couldn't have said it better myself."

Invacare®  
MNS400  
Alternating  
Mattress  
System

# Invacare® MNS400 Alternating Pressure Mattress System



Yes, you can.™